Tortoise[®] Turning & positioning system

Bariatric pressure redistribution system for turning and repositioning

www.molnlycke.com/tortoise

Mölnlycke®

Tortoise® Turning & positioning system



Product and patient alignment

- Place the Mölnlycke Tortoise Turning and positioning pad (hereafter referred as the “pad”) on the bed with the spine- and pelvis print facing up and with a draw sheet on top. Make sure to align the patient with the center line and ensure that the patient’s sacrum is on top of the pelvis printed on the pad. If head of the bed needs to be elevated, confirm that patient’s sacrum is aligned with the lateral sacrum line that goes across the middle of the pad and the hinge joint on the bed.

Following actions require two or more caregivers to perform.

TUCK!

- On both sides of the patient, locate the white handles underneath the pad nearest to the tail. Use these white handles to lift the patient slightly, creating a slight bend at the patient’s knee. Next, take the red handles and tuck the tail under the pad towards the patient’s sacrum. The black non-skid material on the red striped tail should be tucked under the pad, against the bed fitted sheet. No red prints or red handles should be visible after engaging the tail.
- Tuck the black flaps under the pad itself to create an envelopment effect beneficial for pressure redistribution. By tucking the black flaps under the pad, air is displaced creating greater envelopment. This also prevents a sliding motion of the patient. The black flaps can also be tucked under the bed mattress if needed.
- Raise side rails to reduce risk of falls.

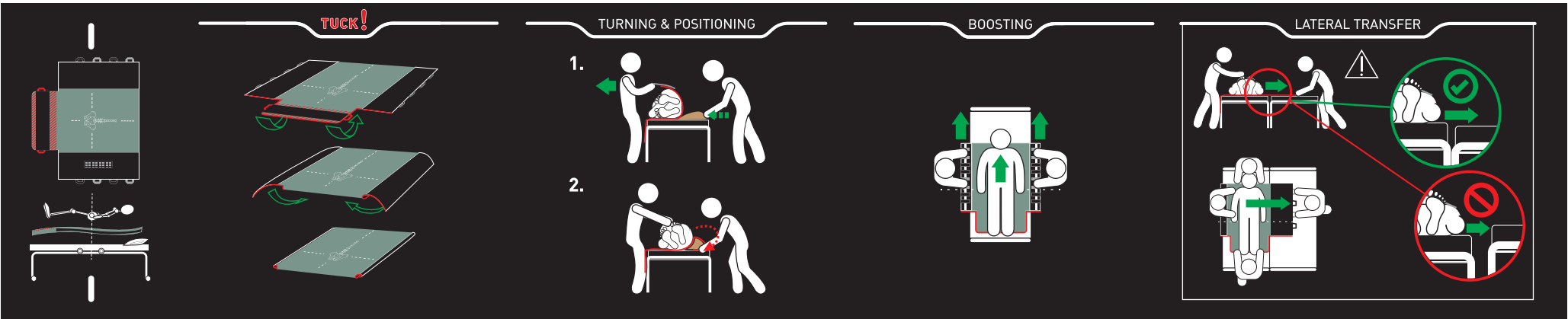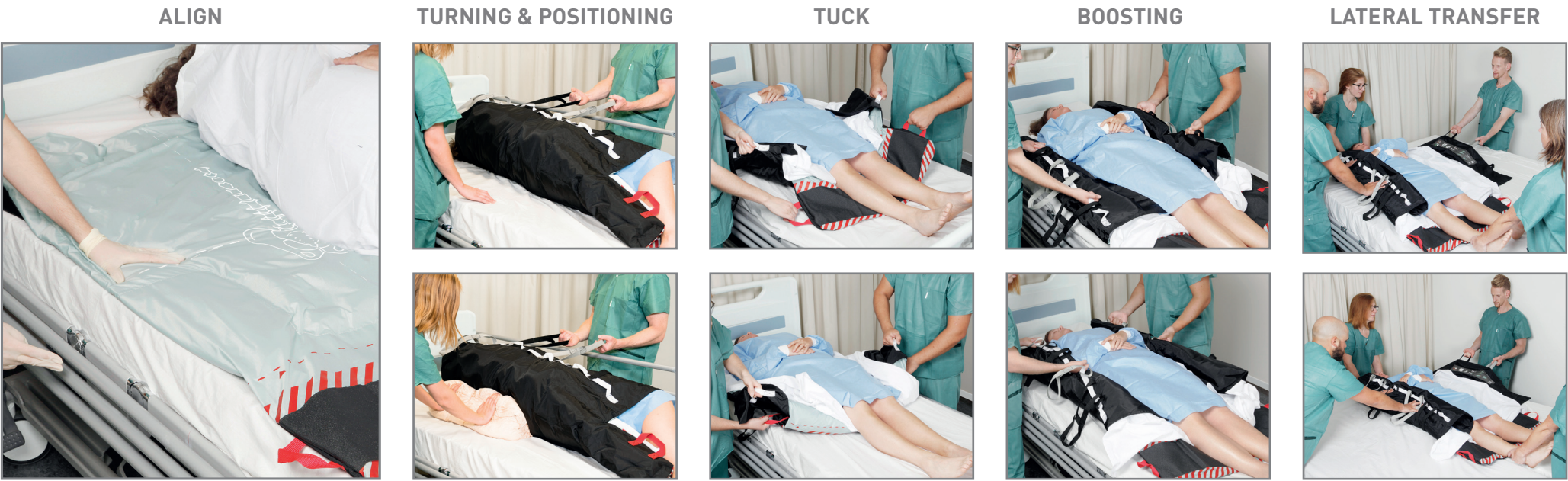TURNING & POSITIONING

- Let one or more caregiver grasp two or more of the long handles on opposite side, palms facing upwards shoulder width apart. To turn the patient, slowly pull the pad towards you until appropriate angle is achieved. Let a second caregiver place the Mölnlycke Z-Flo Fluidized positioner between the pad and the mattress. Mold the Z-Flo towards the body from the scapula down towards the buttock. The Z-Flo can be molded further towards the body to increase the degree of turn.
- Allow the patient to ease back onto the pad. This displaces the air inside the pad and helps to redistribute pressure on the lower hip. Tuck the black flap around the Z-Flo to secure and lock it in place. Tuck the opposite black flap under the pad to create the envelopment effect or tuck it under the bed mattress to prevent a sliding motion of the patient.

BOOSTING

- Locate the white boosting handles on both sides under the black flaps. Ensure each caregiver are using corresponding handles on each side for an effective boost.
- With arms at hip distance and weight evenly distributed, use appropriate ergonomics to boost the patient up in the bed. If appropriate, tilt the bed in Trendelenburg position to facilitate the boost.
- Follow institutional protocol for head, feet and arm positioning. Use the white and the red handles to engage the tail by tucking it under the pad towards the sacrum and tuck the black flaps under the pad (see the TUCK! section above). The black non-skid material should be tucked under the pad, against the bed fitted sheet. No red prints or red handles should be visible after engaging the tail.

LATERAL TRANSFER

The Mölnlycke Tortoise Turning and positioning pad may be used as a lateral patient transfer device. It is not recommended for transfer of patients with documented or suspected instable spinal cord injury. Remember to always follow your facility’s safe patient handling policies and procedures when repositioning patients. Make sure the bed brakes are locked on both beds. Ensure the surface the patient is being transferred to is at same height or slightly lower than the surface they are being transferred from.

- Untuck the tail for an easier transfer.
- If desired, position a lateral sliding aid between the original surface and the end point surface.
- Position members of the patient transfer team on both sides of the patient to prevent the beds from moving apart. Position one caregiver at the patient’s head to support the head and neck during transfer and one caregiver at the patient’s feet to avoid drag if needed.
- The pulling caregiver(s) should grasp two or more of the long handles located on the side of the pad and use slow pulling movement to begin to slide the patient to the end point surface. A heavier patient may require additional pulling caregivers distributed over the four long handles. The white handles alone should never be used for lateral transfer.
- Simultaneously coordinate to gradually help push the patient towards the destination surface while the remaining caregivers are using slow movement to steadily pull the patient toward the destination surface.

en Mölnlycke® Tortoise® Turning and Positioning System

Bariatric pressure redistribution system for turning and repositioning

Always follow your facility’s safe patient handling policies and procedures when repositioning patients.

When using the Mölnlycke Tortoise Turning and Positioning System all actions require two or more caregivers. Use as directed by appropriate clinician and please read all of the instructions before use.

Product Description

Mölnlycke® Tortoise® Turning and Positioning System consists of a positioning pad with a low-pressure air chamber, which adapts to the patient shape by displacing air as the patient is repositioned, providing sacral offloading. It also consists of a rip-stop cover with reinforced handles, allowing caregivers to easily turn and reposition the patient while maintaining proper body mechanics. The system includes a Mölnlycke® Z-Flo™ Fluidized Positioners, assisting the positioning pad in sacral offloading by maintaining neutral body alignment and desired therapeutic angle of patient. The Mölnlycke Tortoise Turning and Positioning System is designed to remain under the patient for continuous support during the length of care.

Intended Use

Mölnlycke® Tortoise® Turning and Positioning System is designed to provide patient turning and repositioning, which may be used as part of a prophylactic therapy to help prevent pressure ulcers.

Indications for use

- Patient turning and repositioning
- Sacral offloading
- Lateral transfer

Contraindications

The Tortoise Turning and Positioning System should not be used with patients who have thoracic, cervical or lumbar fractures.

Intended patient population

Intended for Adults, Teenagers and Pediatrics in hospital patient population.

Lifetime

The product can be used with one patient until the product has deteriorated and does no longer fulfil its intended use. The product expiry date is printed on the packaging labels.

Weight limits

Max 400 kg (880 lbs)

Imaging

The Mölnlycke Tortoise Turning and Positioning System (including the Mölnlycke® Z-Flo™ Fluidized Positioner) is MR Safe and can also be used during CT. The pad can be left under the patient during X-ray, but the Z-Flo itself should not be present in the imaging field.

Product features

Mölnlycke Tortoise Turning and Positioning System consists of:

- Mölnlycke Tortoise Turning and Positioning pad
 - A red striped tail with black non-skid material
 - Black and grey handles for turning & lateral transfer
 - White handles for boosting
- Mölnlycke® Z-Flo™ Fluidized Positioner for patient positioning

Instructions for use

Product and patient alignment

- Place the Mölnlycke Tortoise Turning and Positioning pad (hereafter referred as the “pad”) on the bed with the spine- and pelvis print facing up and with a draw sheet on top. Make sure to align the patient with the center line and ensure that the patient’s sacrum is on top of the pelvis printed on the pad. If head of the bed needs to be elevated, confirm that patient’s sacrum is aligned with the lateral sacrum line that goes across the middle of the pad and the hinge joint on the bed.

Following actions require two or more caregivers to perform. TUCK!

- On both sides of the patient, locate the white handles underneath the pad nearest to the tail. Use these white handles to lift the patient slightly, creating a slight bend at the patient’s knee. Next, take the red handles and tuck the tail under the pad towards the patient’s sacrum. The black non-skid material on the red striped tail should be tucked under the pad, against the bed fitted sheet. No red prints or red handles should be visible after engaging the tail.
- Tuck the black flaps under the pad itself to create an envelopment effect beneficial for pressure redistribution. By tucking the black flaps under the pad, air is displaced creating greater envelopment. This also prevents a sliding motion of the patient. The black flaps can also be tucked under the bed mattress if needed.
- Raise side rails to reduce risk of falls.

TURNING & POSITIONING

- Let one or more caregiver grasp two or more of the long handles on opposite side, palms facing upwards shoulder width apart. To turn the patient, slowly pull the pad towards you until appropriate angle is achieved. Let a second caregiver place the Mölnlycke Z-Flo Fluidized Positioner between the pad and the mattress. Mold the Z-Flo towards the body from the scapula down towards the buttock. The Z-Flo can be molded further towards the body to increase the degree of turn.
- Allow the patient to ease back onto the pad. This displaces the air inside the pad and helps to redistribute pressure on the lower hip. Tuck the black flap around the Z-Flo to secure and lock it in place. Tuck the opposite black flap under the pad to create the envelopment effect or tuck it under the bed mattress to prevent a sliding motion of the patient.

BOOSTING

- Locate the white boosting handles on both sides under the black flaps. Ensure each caregiver are using corresponding handles on each side for an effective boost.
- With arms at hip distance and weight evenly distributed, use appropriate ergonomics to boost the patient up in the bed. If appropriate, tilt the bed in Trendelenburg position to facilitate the boost.
- Follow institutional protocol for head, feet and arm positioning. Use the white and the red handles to engage the tail by tucking it under the pad towards the sacrum and tuck the black flaps under the pad (see the TUCK! section above). The black non-skid material should be tucked under the pad, against the bed fitted sheet. No red prints or red handles should be visible after engaging the tail.

LATERAL TRANSFER

The Mölnlycke Tortoise Turning and Positioning pad may be used as a lateral patient transfer device. It is not recommended for transfer of patients with documented or suspected instable spinal cord injury. Remember to always follow your facility’s safe patient handling policies and procedures when repositioning patients. Make sure the bed brakes are locked on both beds. Ensure the surface the patient is being transferred to is at same height or slightly lower than the surface they are being transferred from.

- Untuck the tail for an easier transfer.
- If desired, position a lateral sliding aid between the original surface and the end point surface.
- Position members of the patient transfer team on both sides of the patient to prevent the beds from moving apart. Position one caregiver at the patient’s head to support the head and neck during transfer and one caregiver at the patient’s feet to avoid drag if needed.

- The pulling caregiver(s) should grasp two or more of the long handles located on the side of the pad and use slow pulling movement to begin to slide the patient to the end point surface. A heavier patient may require additional pulling caregivers distributed over the four long handles. The white handles alone should never be used for lateral transfer.
- Simultaneously coordinate to gradually help push the patient towards the destination surface while the remaining caregivers are using slow movement to steadily pull the patient toward the destination surface.

Items to consider

- SACRUM:** The Patient’s sacrum should always be aligned with the printed pelvis.
- EDGES:** Black edges should always be at the bed edge.
- TAIL:** The tail should always be tucked and no red prints or red handles should be visible.

Cleaning Instructions

Use hospital approved disinfectant or household antimicrobial cleaner / disinfectant according to the label instructions. Do not use bleach. Extended use of alcohol- based disinfectants may reduce the life of the product.

Warnings

- Do not use on patients one year of age or younger due to risks of suffocation.
- Do not puncture or use sharp objects near system.
- Do not use product with lifts.

Precautions

- Do not X-ray through the adult Z-Flo.
- Keep the handles out of imaging area during X-ray
- Handles should be tucked away during CT and MRI to avoid entangling.
- The white handles alone should never be used for lateral transfer.
- No maintenance required. Not user repairable.
- Replace promptly if punctured or damaged.
- If the Z-Flo is punctured, place it in a sealed plastic bag to avoid spillage of the fluidized material.
- Disposal should be handled according to local environmental procedures.
- The use of a linen is required between the patient’s skin and the pad.

Other Information

- The device is for single patient use only.
- Product provided non-sterile. Do not sterilize.

If any serious incident has occurred in relation to the use of Mölnlycke® Tortoise® Turning and Positioning System it should be reported to Mölnlycke Health Care and to your local competent authority.

Special storage conditions/handling conditions

Mölnlycke Z-Flo Fluidized Positioners should be stored below 40°C (104°F).

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